


Stone Extraction Balloons

Single Use

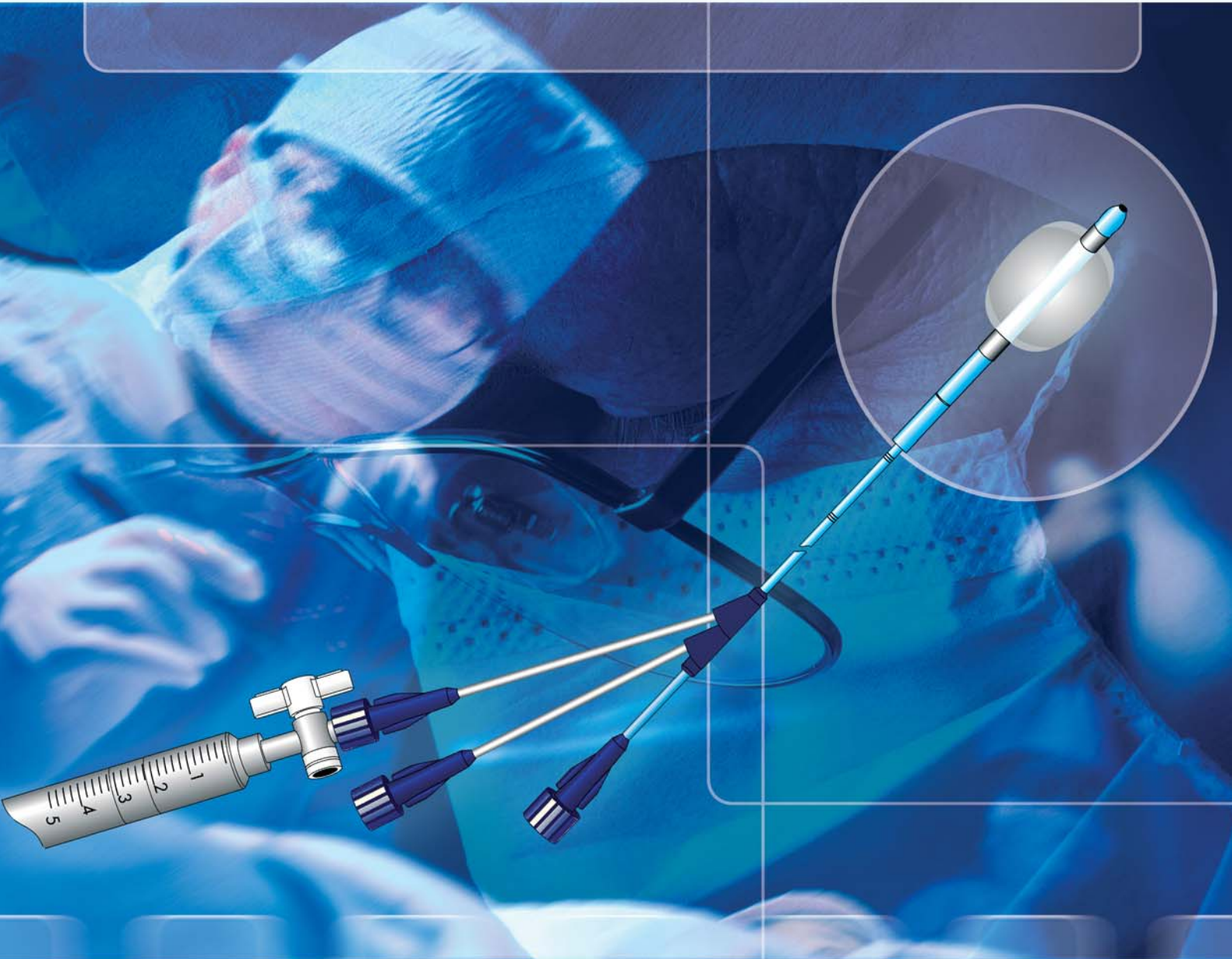

Stone Extraction Balloons

Single Use

Features:

- Double Lumen -----> Guide Wire and Contrast Compatibility
- Three Separate Ports for Balloon Inflation, Contrast Medium and Guide Wire -----> Easy and Effective Handling
- Two-Step Balloon Inflation -----> Two Defined Balloon Diameters per Size
- Metal Markings at the Distal and Proximal End of the Balloon -----> Exact Positioning
- One Centimeter Markings (30 cm) -----> Macroscopically Controlled Penetration Depth
- Removable Stylet (1.75 m) -----> Kink Free Introduction into the Working Channel
- Single Use -----> Safety and Convenience

Benefits:

Caution: This Product Contains Natural Rubber Latex Which May Cause Allergic Reactions.

Order information:

Catalog No.	French size	Length	∅ Balloon	Max. Air Inflation	Max. Liquid Inflation	Accepts Guide Wire
GEB-01-17-025	5.0	200 cm	9-12 mm	3.0 ml	1.5 ml	.021" (0.53 mm)
GEB-01-23-035	7.0	200 cm	14-16 mm	4.5 ml	2.5 ml	.035" (0.89 mm)

Manufactured by:



Medi-Globe GmbH
 Medi-Globe-Strasse 1-5
 83101 Achenmühle, Germany
 Phone: ++49 8032 973 379
 Fax: ++49 8032 973 339
 e-mail: sales@medi-globe.de
 www.medi-globe.de

Distributed by: